



CACI NON SURGICAL
FACIAL TONING

CACI CLASSIC



MICROCURRENT



HYDRATONE

CACI CLASSIC

Non Surgical Facial Toning System

The result of over 20 years research, the award winning CACI facial toning treatment has been voted* the most effective anti-ageing treatment available. The treatment results are often so dramatic that CACI has become known as the "Non Surgical Facelift".

CACI systems are the most advanced and sought after microcurrent machines available. The Classic is the latest generation of the original CACI Mark One facial toning system.

CACI's advanced microcurrent facial toning technology is also featured in our multifunctional Quantum, Ultra and Ultimate systems.

MICROCURRENT

Originally developed for treating facial palsy, microcurrent therapy is now widely used as an aesthetic anti-ageing treatment.

True microcurrent uses a current with an intensity of less than one millionth of an amp (<1mA) that works in harmony with the body's own bio-electrical field. A higher intensity current will cause contractions of the facial muscles which can actually worsen lines and wrinkles, whereas microcurrent stimulation does not contract muscles, but instead works by a process called muscle re-education to achieve results.

In the 1960's Drs Baker & Doyle demonstrated that effective long-term muscle re-education could only be achieved using frequencies less than 1 Hertz. All CACI microcurrent systems use these ultra low frequencies to re-educate facial muscles.



Non Surgical
Treatments To Tone
The Face & Soften
Lines & Wrinkles



MICROCURRENT BENEFITS

- Lifts and tones facial muscles
- Reduces the appearance of lines and wrinkles
- Increases blood and lymph circulation
- Enhances the penetration of skin care formulations
- Increases the production of collagen and elastin
- Increases protein synthesis and cell membrane transport
- Increases mitochondria activity and adenosine triphosphate (ATP) production



FACIAL TONING

CACI gently helps tighten and tone sagging muscles whilst smoothing the appearance of fine lines and wrinkles. The dual-tipped probes increase the effect on the muscle and when combined with the T-bar have dramatic firming results on the jaw line.


An impressive feature of CACI is the simultaneous use of pads and probes, which enable the face and neck areas to be treated at the same time.



HYDRATONE

Hydratone combines active microcurrent rollers with Hydro Mask, a unique, electrically conductive silicon gel mask that has been infused with powerful hydrating properties.

The rollers gently massage the face causing the mask to become electrically charged so that the whole face is bathed with rejuvenating energy. The current delivered by the rollers recharges and replenishes the bio-electrical energy within the facial muscles, restoring tone and firmness.

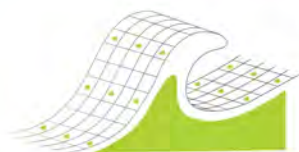

THE CACI DIFFERENCE

All CACI systems deliver microcurrent energy using ultra low frequencies and the power of a unique digitally simulated tidal wave.

This special digital waveform optimises results and comfort by gradually increasing in intensity to overcome skin resistance (Ohm's Law). It then rises to a peak in order to deliver the full power of the microcurrent energy deep into the muscle.


INCREASE YOUR REVENUE WITH CACI CLASSIC

When you buy the CACI Classic, you buy more than just a treatment system, you become part of a winning team.


Share in the success of the CACI brand name and increase your revenue.

*Based on a 6 day week and 1hr treatment at £40

WEEKLY REVENUE

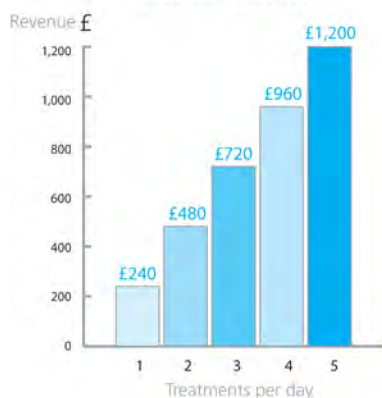

ANNUAL REVENUE

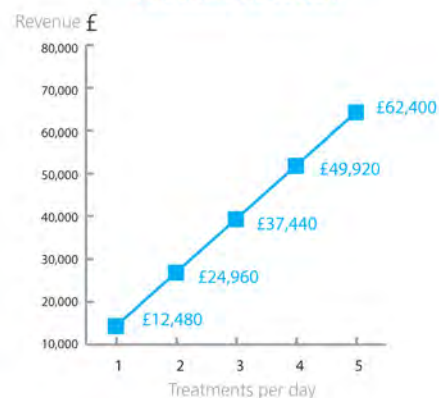



MICROCURRENT



BEFORE



AFTER

During each 1-hour treatment session all 32 facial muscles are gently re-educated.

Microcurrent energy is transmitted to the facial muscles to restore tone and elasticity and to reduce the appearance of lines and wrinkles.



53 PRESET PROGRAMS

TREATMENTS

- Facial Toning ✓
- Lines / Wrinkles ✓
- Eye Bags ✓
- Skin Hydration ✓
- Neck Treatments ✓
- Cellulite (with ECM)* ✓
- Hand Treatments ✓

* Optional Attachment

ELECTRO BUDS



CACI systems utilise specially developed serum filled cotton buds to enhance treatment results. The serum contains Hexapeptides and powerful anti-ageing ingredients to provide an instant face firming effect. Electro buds are also perfect for retail.

ELECTRO CELLULITE MASSAGER (ECM)



The ECM is an optional attachment that can be used in conjunction with all CACI systems that incorporate microcurrent technology. The ECM attachment combines the therapeutic benefits of massage with the effectiveness of microcurrent therapy.

Conductive rollers are electrically charged with microcurrent to gently knead tired shoulders, soothe aching limbs and firm sagging muscles. The stimulatory action of the ECM tones muscles and reduces inches by improving lymph drainage and breaking down fatty deposits.

“ The training and the launch of the CACI machine at the salon was second to none and we ended up selling 7 courses. ”

Laura Ludlow -
Nicetouch,
Chipping Norton





Touch Screen Interface

TECHNICAL SPECIFICATION

Dimensions(cm) 45W x 42D x 15H
Touch Screen 10.4in
Weight 9Kg

Supply Voltage 60-240 volts
Output Voltage: 30 volts max (loaded)
Consumption 110 vA
Safety Class 1 type BF in accordance with IEC 601 and IEC 601-2-10
Channels: 4 independent

MICROCURRENT
Polarity: positive, negative & alternating
Range: 20 uA to 640 uA
Frequency: 0.1Hz - 800Hz
Waveshapes: 8
Mode: Automatic Sequential

INSIDE THE BOX

ATTACHMENTS

- 1x Pair Micro Probes
- 1x ETR Lead
- 1x Mains Lead
- 1x T-Bar
- 1x Pair Straight Bars
- 1x Pair Hydratone Face Rollers
- 1x Pair Double Tip Probe Inserts

CONSUMABLES

- 3x Self-Adhesive Pads -Small (Set of 4)
- 1x Self-Adhesive Pads -Large (Set of 4)
- 1x Collagen Revital Serum (30g)
- 1x Essential Cleansing Milk (300ml)
- 1x Eye Wrinkle Gel (100g)
- 1x Firming Gel Concentrate (300ml)
- 1x Hydro Moisturising Cream (300ml)
- 1x Revitalising Toning Fluid (300ml)
- 1x Pack Electro-Buds (Per 100)
- 20x Hydro Masks

PROMO/MARKETING MATERIAL

- 1x A1 Poster
- 1x A5 Client Leaflets (Per pack)